

Caring for you

Contacts

Headquarters

Neusoft Medical Systems Co., Ltd.
No.177-1 Chuangxin Road, Hunnan District,
Shenyang, Liaoning, China 110167
Email: zhang-dan@neusoftmedical.com

Asia & Oceania

Neusoft Medical Systems Co., Ltd.
No.177-1 Chuangxin Road, Hunnan District,
Shenyang, Liaoning, China 110167
Email: yanghw@neusoftmedical.com

Africa

Neusoft Medical Systems (Africa) Co., Ltd.
D1,Ground Floor,Morningside Office Park,
Ngong Road, Nairobi, P. O. Box 22288-00505, Kenya
Email: yu.xm@neusoftmedical.com

Europe

Neusoft Medical Europe GmbH
Mergenthaler Allee 45
65760 Eschborn, Germany
Email: morinl@neusoftmedical.com

Middle East

Neusoft Medical (MENA) FZ- LLC
No. 705/706, Building 26, Al-Baker Building
Dubai Healthcare City, P.O.BOX 115321, UAE
Email: liuwanj@neusoftmedical.com


North America

Neusoft Medical Systems, USA Inc.
14425 Torrey Chase Blvd, Suite 100
Houston, TX 77014, USA
Email: christopher.mchan@us.neusoftmedical.com

South America

Neusoft Medical Peru S.A.C.
Calle Los Conquistadores
175a, San Isidro, Lima 27, Peru
Email: liuba@neusoftmedical.com

www.neusoftmedical.com/en

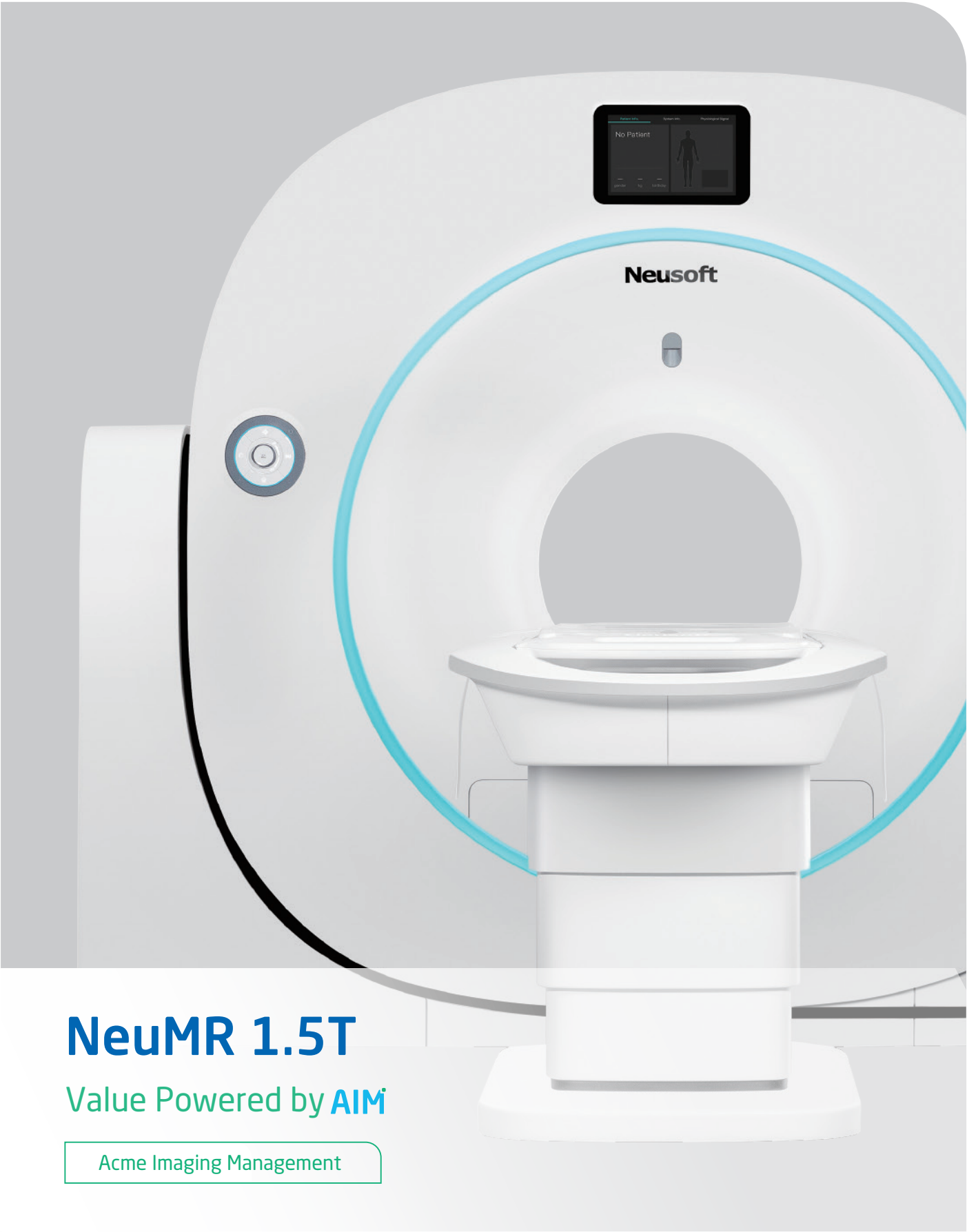

NeuMR 1.5T
Value Powered by AIM

Acme Imaging Management

NeuMR 1.5T/Y1911




NeuMR 1.5T

Value powered by **AIM**, NeuMR 1.5T is designed to offer friendly interactive experience and streamlined workflow while delivering robust and exquisite imaging to customers.

AIM platform is the engine of NeuMR 1.5T which includes newly designed RF coils, table, cover and UI for better user experiences; ARM gradient coil for stable image quality.

NeuMR 1.5T provides comprehensive and quantitative applications, such as Total Neuro, Large-scale DIXON, BrainQuant, HDDWI&DTI, MRS, ASL, etc. Neusoft **AIM** MRI expands your clinical capabilities.

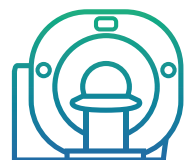



AIM Platform for Streamlined Workflow

AIM coils

- Right and easy connection
- Integrated and combined imaging
- High element density
- Compatible with NeuSENSE
- Integrated with K-fiber



Easy coil handling & fast patient positioning



- Flexible
- Combinable
- High element density
- Lightweight
- General purpose
- Optimized coverage
- Reliable
- RightConnect
- EasyConnect

AIM Torso coil




AIM GP Flex M coil



AIM GP Flex S coil



AIM Platform for Friendly Interactive Experience

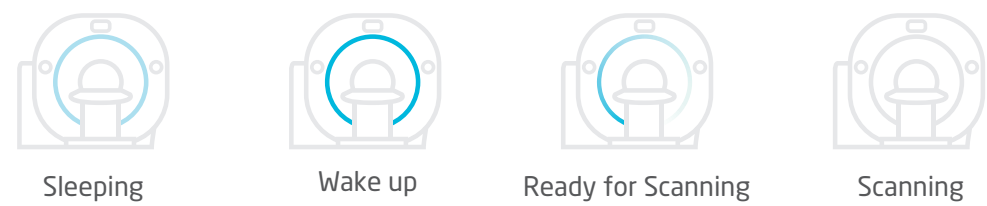

-  Anthropomorphic light belt indicating system status
-  Large LCD displaying system and patient informations
-  Ergonomic control panel for easy operation
-  Integrated **AIM** coils for streamlined workflow
-  Wide table top for comfort and safety
-  Large side cover for service convenience




Aesthetic layout design




Intuitive Designed UI for Friendly Use

Graphical posture for body position confirmation
Graphical anatomy selection for easy protocol setting

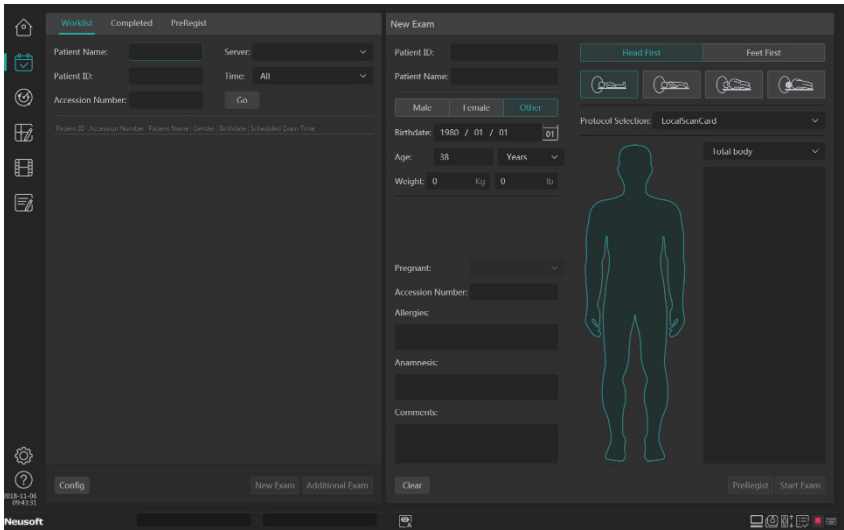

Mosaic image preview for selected information viewing in advance
Shortcut of Data analysis tools for easy and fast operation

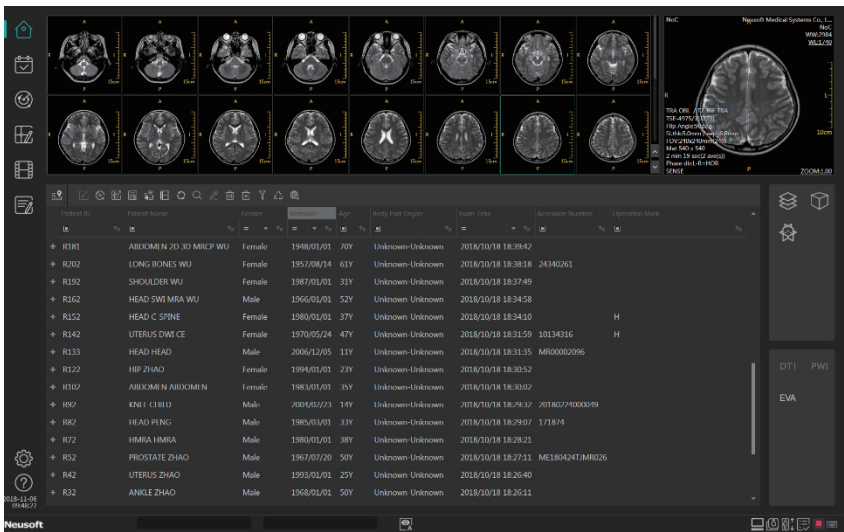


ARM GC for Stable Image Quality

Accurate Temperature Management Gradient Coil

Multi-channel cooling design with multi-point temperature control enables real-time gradient coil temperature control, which realizes high spatial linearity and high duty cycle performance. Using latest eddy current compensation technology, the gradient system delivers fidelity gradient waveform, hence, stable image quality




With ARM GC, NeuMR 1.5T realizes excellent PC MRV and high quality DWI images.

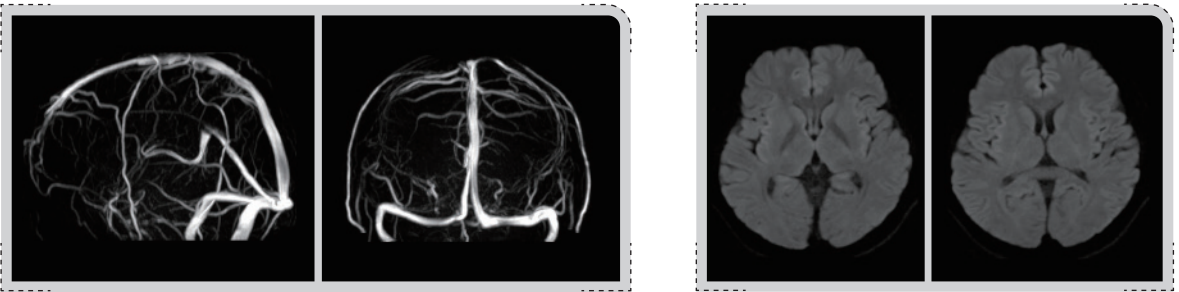

Comprehensive Clinical Applications




Neuro ▼

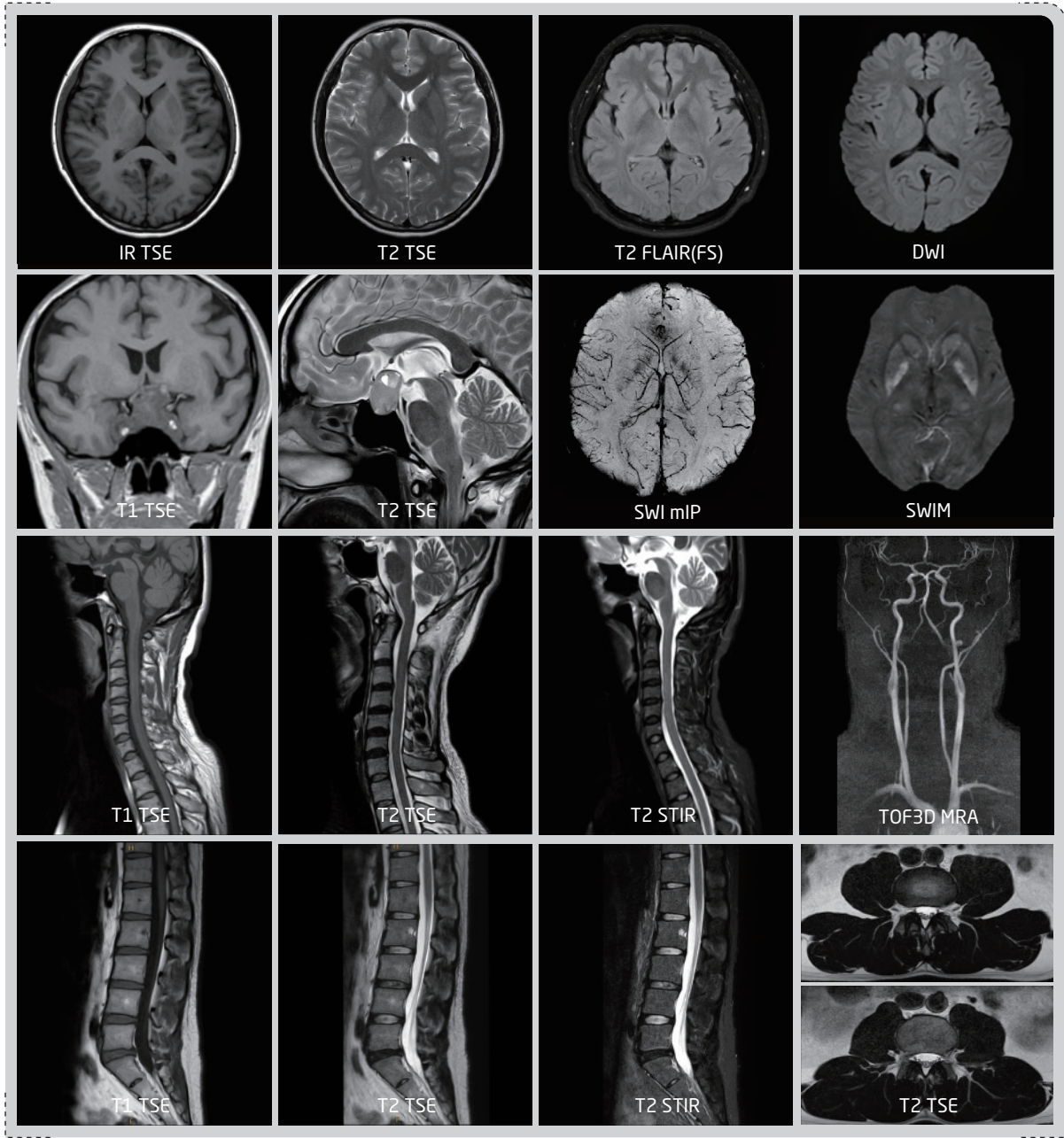


Comprehensive Clinical Applications



Total Neuro ▼

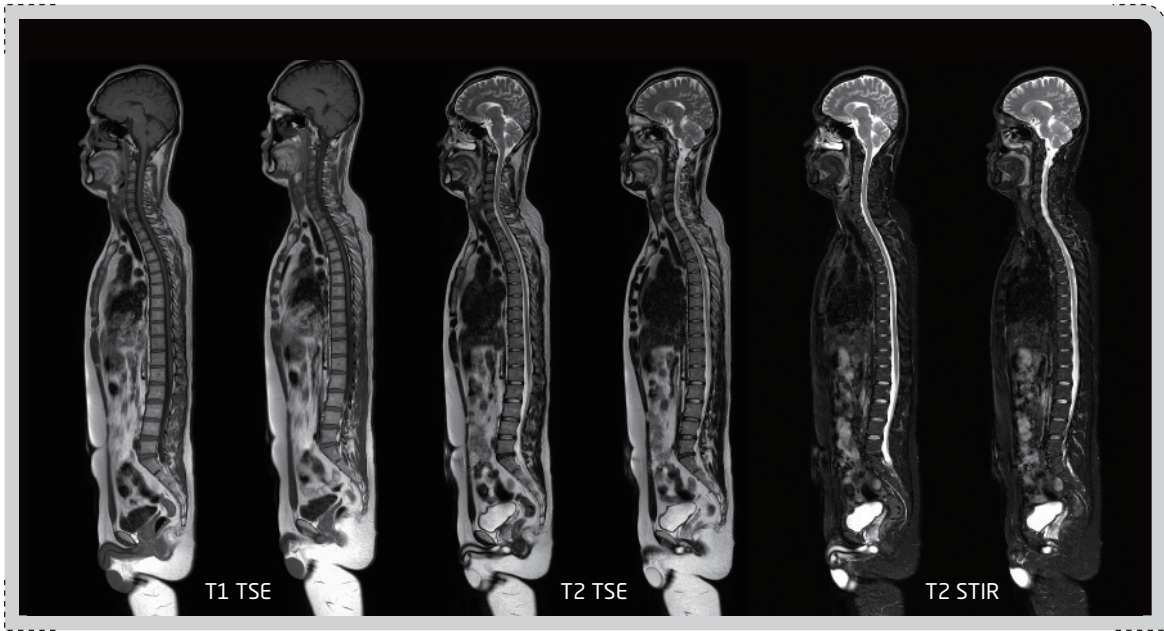


Large-scale DIXON* ▼

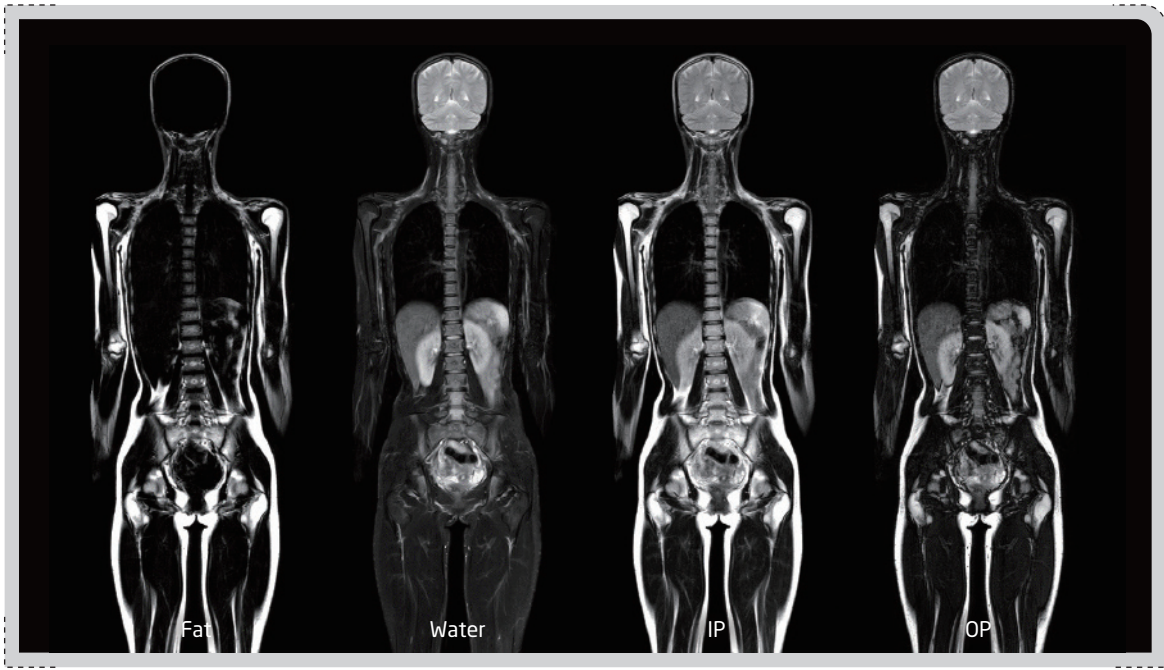

*Option

Body ▼

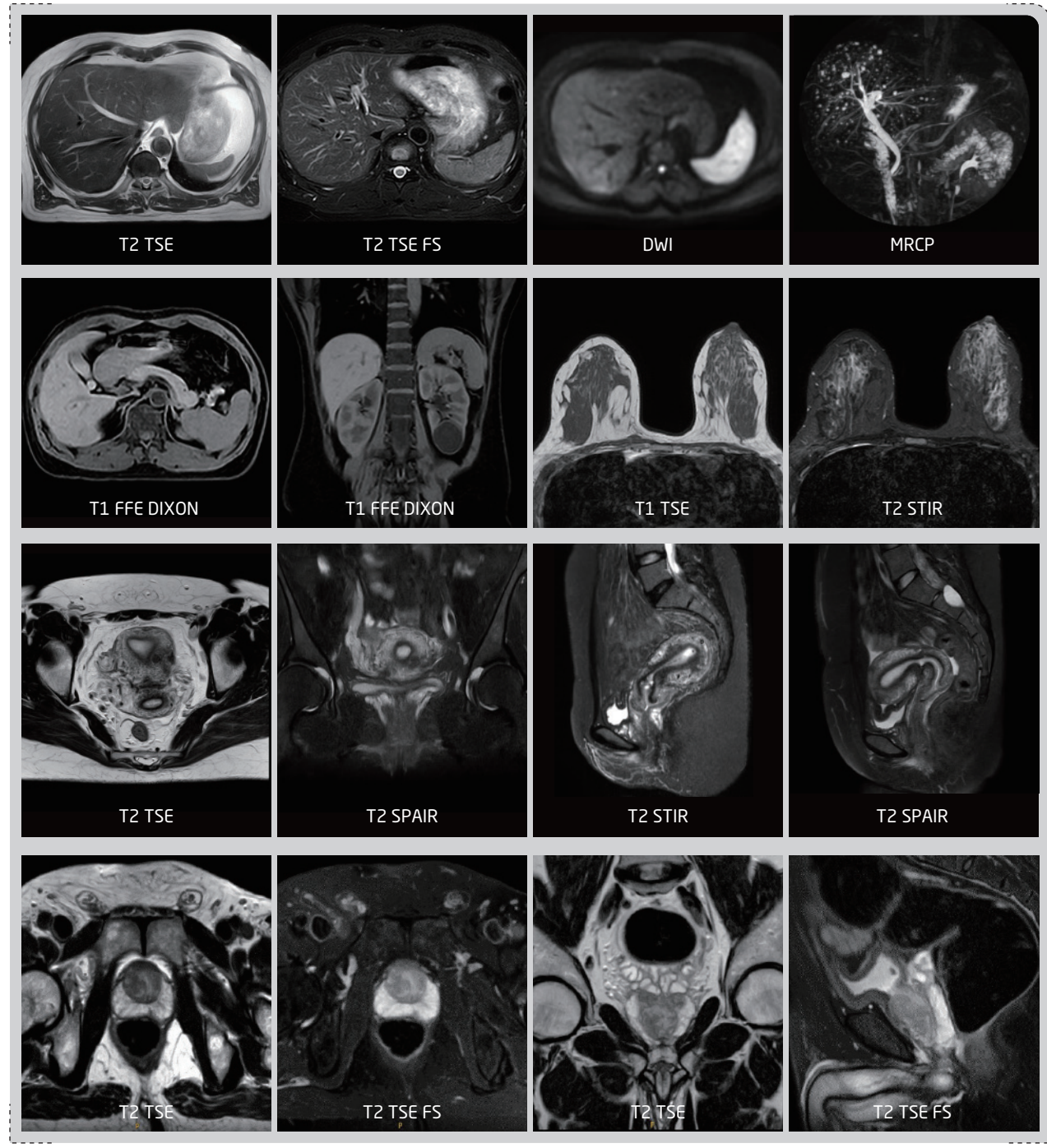


Comprehensive Clinical Applications

MSK ▼





DIXON ▼

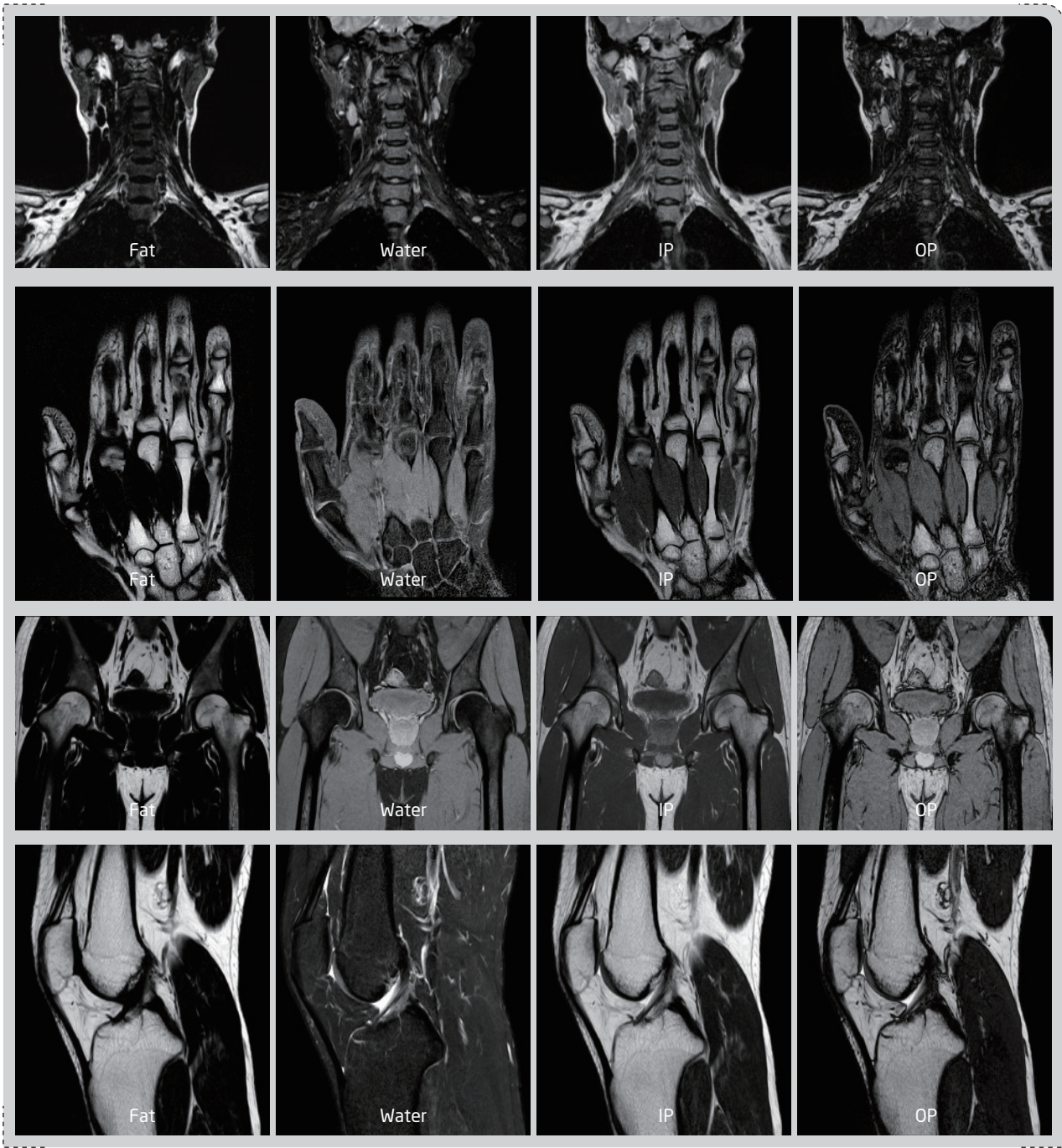


Advanced and Quantitative Applications

BrainQuant ▼

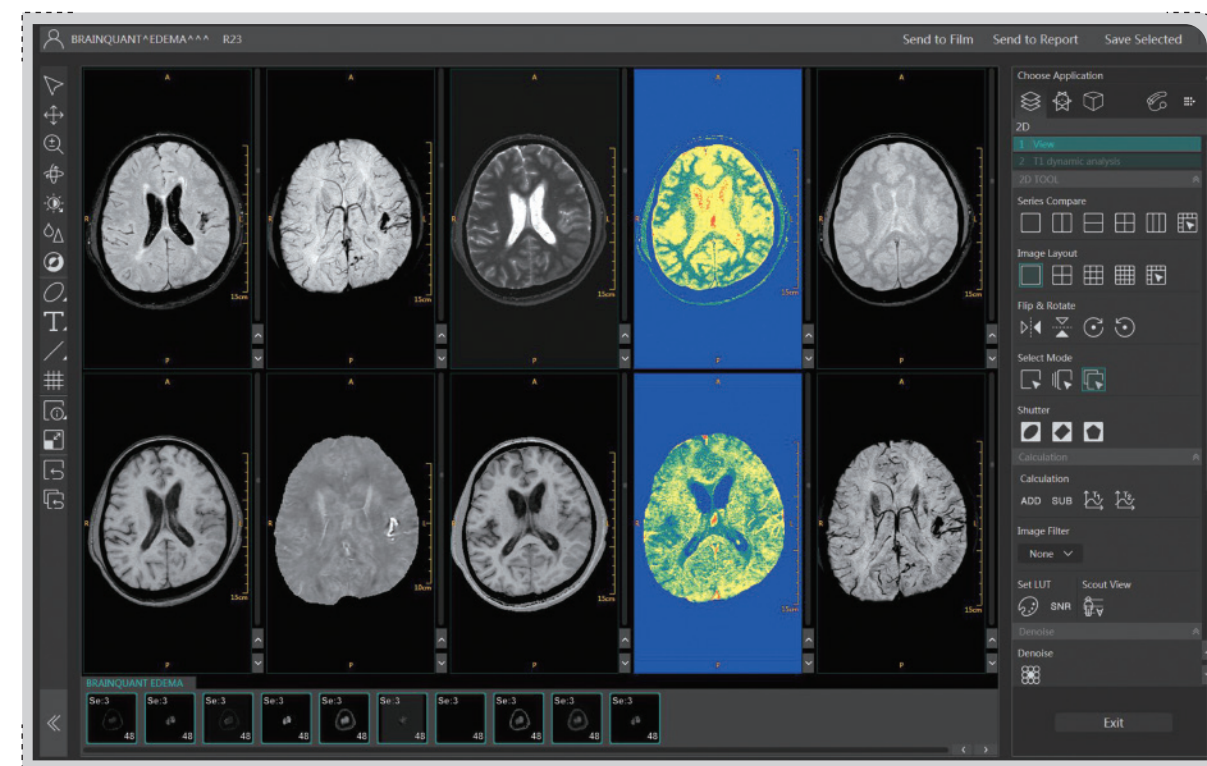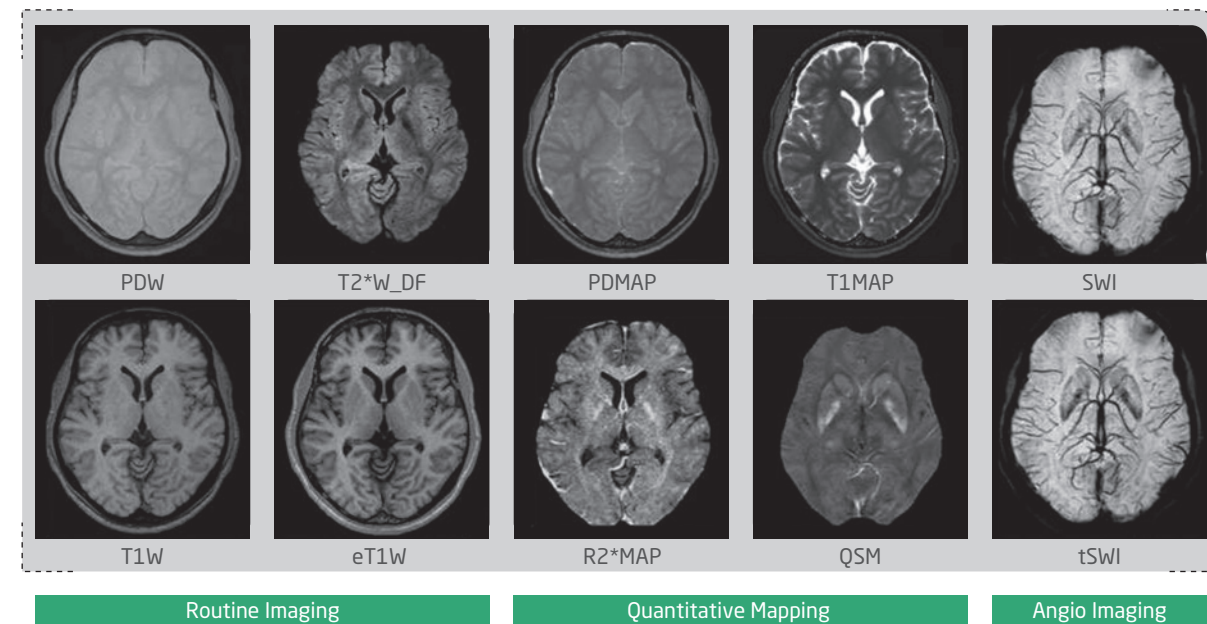


HDDWI & DTI ▼

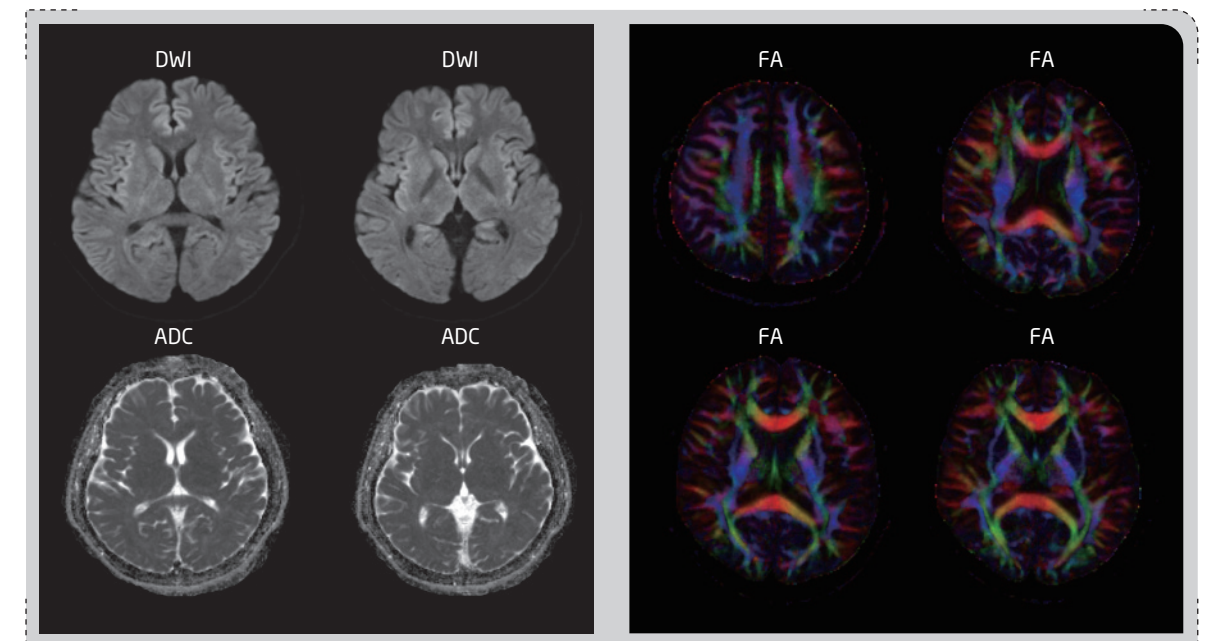

b value=1000
Matrix: 192×192

Direction No. =32
Matrix: 192×192

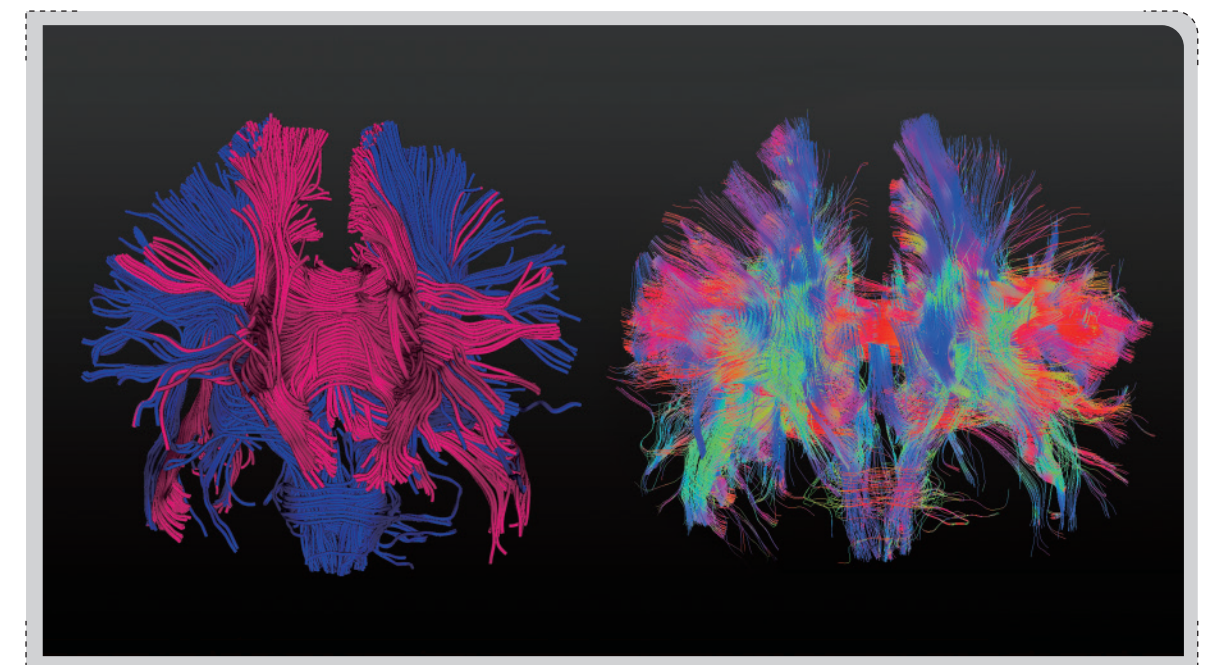



RAPID Service System - Power in Your Hands

Simplifying customer report process
Improving service efficiency
Maximizing system runtime

Neusoft Global Service & Logistics Network

One-click for service report

Online remote service

Predict spare parts lifetime



After-sales service and support

- Remote service capabilities bring Neusoft expertise to you IMMEDIATELY, no matter where you are!
- Identifying and correcting PROMPTLY and PROACTIVELY, minimizing downtime and patient inconvenience.
- Global logistics network enables fast response regarding parts and supplies.

* Note: The contents of this publication and the listed parameters are for reference only and not intended as legal offers or commitments. Neusoft Medical Systems reserves the right to modify the contents, design, specifications and options described herein without prior notice, and will not be liable for any consequences resulting from the use of this publication. Please contact your local Neusoft sales representative for the current information. The specific sales product configuration is subject to the actual contract signed by Neusoft.

* Not available in the United States.

HQ

Subsidiaries

Service centers

Logistics center

Training center

Office